



TopCEM

Dual Cure Resin Cement

Upgraded stability,
consistent versatility.

An all-rounder for all types of
restorations and materials
cementation.



All For Unsurpassed Stability.




Vetted in a series of long-term trials, VinciSmile rolls out PTHP, an innovative redox self-stabilizing technology, to the whole dentistry.

By controlling the oxidation state, PTHP core technology not only maintains the long-term stability of the physical and chemical properties, but significantly reduces the impact of factors such as high temperature and long storage times.

High
Temperature
Resistance

PTHP technology mitigates the high temperature sensitivity of the cement, thereby addressing potential premature curing issues during transportation and enhancing the overall product performance stability.

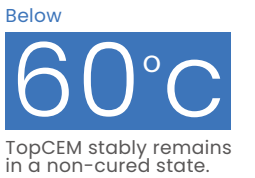


Stable
Handling
Performance

PTHP technology lowers the requirements for transportation conditions. It ensures color consistency before and after curing and provides long-term stability in working and curing times, finally achieving long-lasting and reliable cementation.

Impressive
Color
Consistency

TopCEM provides 7 artistic shades and 3 opacity level for indirect restoration. PTHP technology remarkably enhances the color stability, perfectly addressing the issue of color change before and after curing, and achieving optimal aesthetic results.

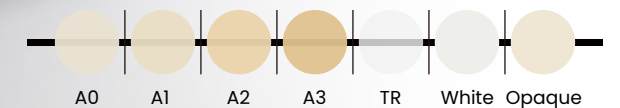

Confidence Drew From Reliability.

TopCEM is exactly what you need for versatile restorations cementation and dental success.

• Excellent Operation Experience

- Sufficient working time ($\geq 150\text{sec}$) delivers high operation flexibility in practice.
- Impressive radiopacity **over 400%** enables easier diagnosis, treatment, and monitoring of dental restorations at recall.

• Optimal Cementation Performance

- A thin film thickness about **13 μm** , ensures the restoration can be fully seated, helps in achieving better marginal integrity, and leads to ideal bonding.
- High bonding strengths are delivered on multiple substrates, decrease the likelihood of fracture or debonding under functional stresses.



• Reliable Longevity

- Ideally low solubility and water absorption significantly help in maintaining TopCEM properties over time, ensuring the durability and longevity of the restoration.
- TopCEM performs out exceptional flexural strength **over 120MPa**, essential for the long-term success and reliability of dental restorations.

| | |
|---|-------------------------------|
| Flexural Strength/ MPa | 123MPa |
| Water Absorption/ $\mu\text{g}/\text{mm}^3$ | 16 $\mu\text{g}/\text{mm}^3$ |
| Solubility/ $\mu\text{g}/\text{mm}^3$ | 0.6 $\mu\text{g}/\text{mm}^3$ |

* Data from VinciSmile Lab.

More Than A Cement.



TopCEM is delivered as a fully aligned and all-round system designed to simplify indirect restoration workflow and elevate your success in dental practice.



Tooth Etching

P-Etchant Phosphoric Acid Etching Gel, with ideal consistency for accurate placement, eliminates the surface dirt and roughens the surface to enhance bonding strength.



Adhesive Application

HugeBond Universal FliPro Light Cure Dental Adhesive, a leading innovation of 8th generation bonding agent, containing MDP monomer, meets all your expectations for a universal bonding.



Restoration Etching

HF-Etchant Hydrofluoric Acid Etching Gel, a reliable etchant of 9% concentration, is used for etching glass ceramics restorations.



Priming

TopCEM Ceramic Primer Ceramic Coupling Agent, a universal primer for conditioning all common restorative materials, significantly improve the bonding strength between the resin cement and the restorations, especially for glass ceramics.



Cementation

TopCEM Dual Cure Resin Cement, an all-rounder for all types of restorations and materials cementation.



Order Information



| TopCEM Dual Cure Resin Cement | |
|-------------------------------|------|
| A0 | 1*8g |
| A1 | 1*8g |
| A2 | 1*8g |
| A3 | 1*8g |
| TR | 1*8g |
| UO | 1*8g |
| WO | 1*8g |



| TopCEM Elite System | |
|-------------------------------|--|
| HugeBond Universal FliPro | 1*6ml |
| TopCEM Dual Cure Resin Cement | 1*8g |
| TopCEM Ceramic Primer | 1*5ml |
| P-Etchant | 1*3ml |
| Accessories | Intraoral mixers Dispensing tips Mixing pad Micro-brushes |



| P-Etchant Phosphoric Acid Etching Gel | |
|---------------------------------------|--------------------------|
| Normal Package | 1*3ml+10*Dispensing tips |
| Bulk Package | 2*5ml |



| HugeBond Universal FliPro Light Cure Dental Adhesive | |
|--|-------|
| | 1*6ml |



| HF-Etchant Hydrofluoric Acid Etching Gel | |
|--|----------------------------|
| | 1*1.2ml+10*Dispensing tips |



| TopCEM Ceramic Primer Ceramic Coupling Agent | |
|--|-------|
| | 1*5ml |



About VinciSmile®

— Your Partner In Prosthetic Dentistry —

VinciSmile Group LLC is a worldwide dental supplier headquartered in Los Angeles, California USA. We offer a wide range of dental lab and clinical solutions, including NOBIL-DENT Resin Teeth, NOBILCAM CAD/CAM Materials, NOBILTRAY Lab Materials, PERFIT Elastomeric Impression Materials, Clinical Restorative Materials, and Invisible Orthodontic Appliance, etc.

Our company is committed to providing exceptional support services, offering comprehensive counseling and training for employees, agents, and customers. As an FDA-cleared facility, we adhere to strict production quality assurance standards and maintain ISO 13485 certification.



VINCISMILE GROUP LLC

Add: 20524 Carrey Rd. Walnut, CA 91789, USA

Email: marketing@vincismile.com Tel: +1 626 283 5808

Web: www.vincismile.com Shop with us online at store.vincismile.com

   
Follow @VinciSmile